Clever

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

Thank you for your purchase and thank you for the trust you have shown. Read carefully before using the product and keep the instructions for possible future use. In case you give the device to another user, please, also give this user instruction manual.

INFORMATION FOR YOUR SAFETY

1. The power supply of your device must match the originally intended specifications.
2. Do not immerse the device in water. It resists only splashes of water.
3. Do not place flammable objects, explosive substances or dangerous objects near your appliance.
4. Use only the accessories and connectors provided. The use of any other type of accessory not intended for this purpose could damage your device irreparably.
5. Use and keep your device in an environment where the temperature is between -10°C to 40°C .
6. Do not disassemble your device and do not try to repair it yourself.
7. This device is not intended for use by persons (including children) whose physical, sensory or mental capabilities are impaired, or persons without experience or knowledge, unless they have a person with them who is responsible for their safety and surveillance or to provide prior instructions concerning the use of the device.
8. Keep away from children.
9. Do not let the device charge unattended. In the event of a short circuit or mechanical damage to the internal Lithium battery, there is a risk of overheating and fire.
10. Do not use your device if it has been damaged.
11. To avoid the risk of electric shock, unplug your device from the main power supply and all devices before cleaning or when not in use.
12. No guarantee can be applied in case of misuse.

13. The manufacturer declines all responsibility for material damage or damage (physical or moral) due to misuse and / or non-compliance with safety instructions.

Operating instructions

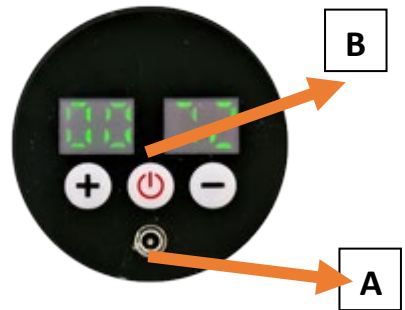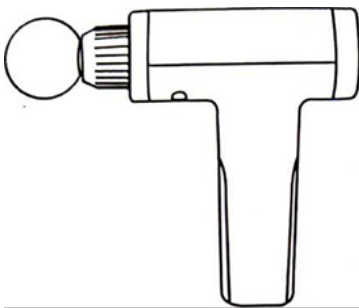1. Connect the charging cable to the socket (A) on the back, then connect the usb to a 5V2A charger
2. The right digital display will light up immediately and show the charge level, when No 99 appears on this display, the device is fully charged and you need to unplug it.
3. Insert the head you want simply by applying it to the hole.
4. Press and hold the power key (B) to turn on the device and the two screens light up, the left screen will be at zero and the right screen will show the battery level.
5. To start the vibrations press the + key, in total you can increase the level by 20 points, each time you press the + will go up one level. Respectively, to decrease the level, press -.
6. To switch off the appliance, hold down the power button for a few seconds. While if you want to pause, press the power key once to enter standby, press again to start.
7. To remove the head simply pull it
8. Always start the program with a low level and gentle movements. 15-30 minutes of massage is enough for each muscle group.

6 massage heads

Ergonomic design

Care for every muscle group

1

Small round head

For joints

2

Long round head

For pain point massage

3

Conical head

For pain point massage

4

Spade head

For pain point massage

5

Spade head

For spine, lumbar spinal
cord massage

6

Round head

For the whole body
muscle massage relaxation

This product is in strict conformity with all the valid directives of European Union applying on this type of electrical appliances.

Protecting the environment

Should you find one day that your appliance needs to be replaced or if it is of no further use to you, think of the protection of the environment:

- 1) Do not dispose your appliance along with the rest of the public waste (this is also the meaning of the shown recycling sign).
- 2) Contact your Public Authorities and they will instruct you of the recycling centers to which your appliances must be disposed.
- 3) Correct disposal of your appliance helps the protection of the environment as well as the recycling of the appliance components.
- 4) The electrical appliances due to their construction materials, if not correctly disposed, may lead to environmental and furthermore health hazards.

LAMD A E-ELECTRONICS P.C.

Peiraios 1 (in stoa Lykourgou), 105 52 Athens
Tel: 210.3220003 e-mail: ecomm@lamda-el.gr

Visit our web site: www.lamda-el.gr